



THE BEGINNING OF THE HERA OF RADIOFREQUENCY



HERA Cellular Coach, the beginning of the new era of radiofrequency. It is the ultimate device for discovering the best tissue appearance. Its way to achieve it is through a bidirectional and respectful dialogue between technology and organism. The treatments reach deep layers to promote tissue regeneration and recovery from cellular biostimulation.

HERA Cellular Coach achieves visible, immediate and lasting results. Its cellular language re-educates and trains the organism.





TECHNOLOGY

HERA Cellular Coach is the only platform that combines **four types of radiofrequency at 448 kHz** and the effectiveness of vacuum therapy. All in a single device to achieve non-invasive, pleasurable and visible results from the first session.



INCREASES VASODILATION
AND OXYGENATION OF TISSUES

STIMULATES
CELLULAR METABOLISM

IMPROVES TISSUE FUNCTION

INCREASES THE PRODUCTION
OF COLLAGEN AND ELASTIN

BALANCES THE **REGENERATIVE PROCESSES**

DRAINS AND DETOXIFIES


PROMOTES LIPOLYSIS

DECREASES FIBROSIS



BIOLOGICAL EFFECTS OF RF

HERA Cellular Coach activates the blood system providing increased oxygenation to the tissue. It contributes to **cellular regeneration and repair.** The synergy of technologies generates the following summative biological benefits on the treated tissues:

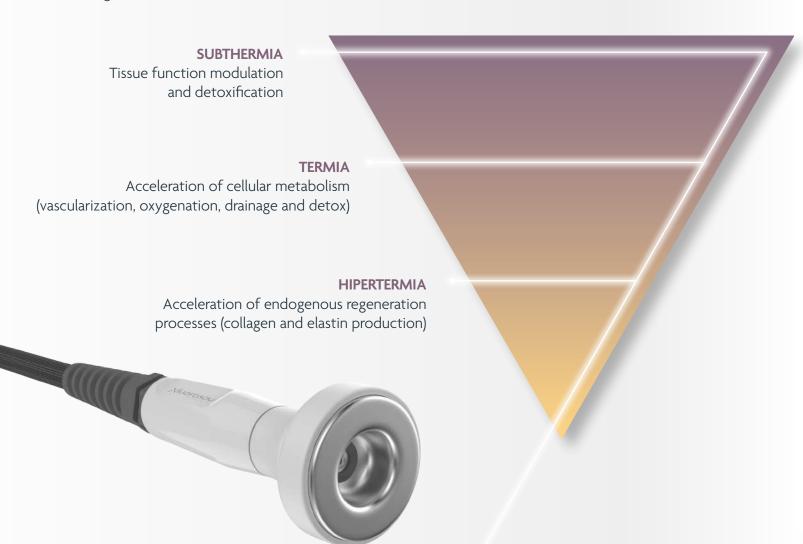







APPLICATIONS AND BENEFITS

FACIAL TREATMENTS:

- Facial rejuvenation
- Skin regeneration
- Wrinkle reduction
- Increased hydration and luminosity
- Increased elasticity

BODY TREATMENTS:

- Localized fat reduction
- Firming treatment for flaccidity
- Treatment for cellulite in all its degrees
- Tired legs

HAIR TREATMENTS:

- Provide vitality and shine
- Stop hair loss
- Revitalize damaged hair
- Densify the hair population

HAND PIECES



HEBEGY hand piece

Combination of vacuum therapy and resistive monopolar radiofrequency. Mobilization, detachment and drainage at its maximum depth.

Ø 30, 55, 75 mm







RESISTIVE MONOPOLAR RF

Indicated for fibrotic and fatty tissue. It allows working the maximum tissue depth. Ø 40, 60, 80 mm



CAPACITIVE RF

Acts on blood circulation regenerating the skin and activating collagen production. Medium depth.

Ø 20, 40, 60, 80 mm



BIPOLAR RF

Highly effective superficial application in the treatment of flaccidity and stretch marks. Ø 30, 70 mm



SPECIAL CURVED ELECTRODES

Designed to facilitate the treatment of certain difficult areas such as double chin, knees, elbows, etc. Thanks to its ergonomic and curved shape. **20, 25, 30 mm**

They are also excellent for capillary areas. The application of radiofrequency generates vasodilation and oxygenation of the scalp. The improvement of irrigation favors the arrival of nutrients and amino acids necessary for the cells of the hair follicles to multiply and regenerate, thus strengthening the hair and favoring its growth.



PROTOCOLS

To facilitate the work of the professional and achieve the best results in the shortest possible time, **Hera Cellular Coach** is based on treatment protocols designed according to the area and the aesthetic alteration to be treated.

FREE PROGRAMME

that allows personalising the treatment to 100% according to the parameters established by the professional.

Treatment
PROTOCOLS
that guarantee
optimal results.

PROGRAMMES COMBINABLE BETWEEN THEM





SMART TIME TECHNOLOGY

by Novasonix

SMART TIME is a unique system that allows the selection of the energy delivery method and temperature for more pleasant treatments, enhancing the cabin experience for even the most sensitive skin.



HERA

TECHNICAL SPECIFICATIONS

GENERAL	Input voltage:	100-240V
	Monopolar output power:	CAP 500VA, RES 200W
	Bipolar output power:	RES 100W
	Frequency:	448 kHz
	Monopolar capacitive electrode:	ø 20, 40, 60, 80 mm
	Monopolar resistive electrode:	ø 40, 60, 80 mm
	Bipolar resistive electrode:	ø 30, 70 mm
	Hebegy electrode:	ø 30, 55, 75 mm
	Display:	10" táctil
	Console weight:	8-9 kg





ADVANTAGES OF HERA CELLULAR COACH

HERA Cellular Coach is the most advanced equipment that easily, intuitively, safely and effectively complements the talent and experience of aesthetics and skin health professionals.



- FREQUENCY 448 KHZ: proven maximum effectiveness in cell and tissue regeneration.
- **CONSTANT EMISSION**: a stable radiofrequency signal that guarantees the development of the natural biological effects on the treated tissues.
- **SMART TIME** temperature control.
- **HEBEGY HANDPIECE**: the only handpiece on the market with vacumm + monopolar RF with interchangeable electrodes to adapt to any area of the body.
- THE ERGONOMIC DESIGN of the new handpieces adapts to the practitioner's hands and different gripping to the practitioner's hands and various handholds with maximum comfort.
- THE DIFFERENT SIZES OF THE ELECTRODES make it possible to work on any area, providing comfort during treatment, even in areas such as knees, elbows and masseter.
- **ELEGANT AND COMPACT:** its avant-garde design results in an elegant and light platform, easy to move and handle.
- **HIGH VERSATILITY:** Hera makes possible a wide range of facial, body and hair treatments on the same platform.
- **HIGH PROFITABILITY:** the return on investment is perceived in the short term when starting treatments.





HEBEGY HANDPIECE

- MONOPOLAR OUTPUT POWER RES 200W
- MAXIMUM ASPIRATION 300 MBAR
- MAXIMUM AIRFLOW 7 LPM
- ELECTRODES 30, 55, 75 MM



Novasonix

There is always a new beauty. You just have to know how to discover it.

We are passionate about beauty and, we believe in it. We know that the skin is the way to achieve it. The gateway to specific care from the inside to reach the maximum splendour on the outside. To this end, we place all our expertise at the disposal of a single goal: to treat each body in a unique way to achieve all the health and beauty it deserves.

All our devices have cutting-edge, non-invasive and safe technologies to satisfy the most demanding beauty professionals and their customers. Proof of this is the constant investment in R&D&I, the inexhaustible desire to innovate and the creation of our patents such as the NSX method.

Our trajectory manufacturing medical-aesthetic equipment has led us to be present in more than 20 countries, basing our development on trust, results, professionalism and the endorsement of the highest European quality and manufacturing certifications.

